INTENDED USE:

The Foldable Transfer Wheelchair “I-TRAVEL” is designed to transport people with mobility issues or limited mobility always with the help of a caregiver. The lightweight aluminum frame and portable design, featuring an ergonomic carrying handle and easy folding mechanism, ensure comfortable transport and use. Finally, the safety belt enhances the user’s safety.

INDICATIONS:

The wheelchair is intended for use both indoors and outdoors while its maintenance, transport, and supervision are handled by a caregiver.

CONTRAINDICATIONS:

The wheelchair should not be used if the user has problems related to balance or perception disorders. The maximum weight capacity, as indicated by the manufacturer, must not be exceeded.

PATIENT GROUP:

- Patients with neurological conditions.
- Patients with mental health conditions.
- Patients with musculoskeletal problems.
- Patients with paraplegia.
- Postoperative patients.
- Patients with chronic cardiorespiratory problems.
- Patients in the rehabilitation phase following lower extremity or pelvic fractures.

INTENDED USERS:

- Individuals who require assistance when transferring to and from the bed, toilet, or another seating surface.
- Wheelchair users, when an alternative means of transportation is required to move them indoors or outdoors.
- Caregivers, healthcare professionals, or other appropriately trained individuals who assist the user with transferring, positioning and using of the wheelchair.
- Elderly individuals

GENERAL SAFETY WARNINGS:

- DO NOT use this product or any available optional equipment until you have first read and fully understood this user manual.
- If you are unable to understand the warnings, precautions, or instructions, contact a healthcare professional, representative, or technical staff before attempting to use the device; otherwise, injury or damage may result.
- Keep the product away from heat sources.
- Be careful when children are nearby, and DO NOT allow children to play with the product.
- Do not exceed the maximum safe load.
- Always exercise extreme caution around moving parts that could cause limbs to become trapped and result in injury.
- Do not attempt to lift the wheelchair by any of its assembled parts.
- Getting in and out of the wheelchair without assistance requires sufficient upper-body strength and mobility.
- Be extra careful with wet and slippery surfaces.

⚠ CAUTION:

- In case of damage or malfunction of your product DOT NOT use the product and please contact the authorized and trained representative.
- If you don’t read this user manual, it is preferable not to use this product or another available part.
- Do not use the product for a purpose not indicated in this manual.
- MOBIAK S.A declines all responsibilities for any consequences resulting from an incorrect use of this product and from unauthorized alteration to the frame of the product.
- MOBIAK S.A reserves the right to change the information contained in this document without previous notice.
- There is a ± 3%. Deviation for the Dimensions of the Products.

BEFORE USE:

- Check if the frame is damaged in order to guarantee a safe use of the product. (There are no cracks or fractures in the frame).
- In case of damage do not use the product and contact your dealer for further instructions.
- Check, if the product is properly assembled and all screws are secure and well screwed on.
- Always check the condition of wear for the parts to ensure the absolute safe use of the product.
- We do not recommend that the product is issued to anyone without first carrying out a full assessment of their specific needs. It is important that the user is given the necessary support by the product as prescribed by a qualified professional.

ENVIRONMENTAL PROTECTION:

If one day you find that your product needs to be replaced or it is no longer working for you, consider protecting the environment:

1) Do not dispose your product along with the rest of the public waste (this is also the meaning of the shown recycling sign).
2) Contact your Public Authorities and they will instruct you of the Recycling centers to which your product must be disposed.
3) Correct disposal of your product helps the protection of the environment as well as the recycling of the product’s components.

INCIDENT REPORTING:

Any serious incident that occurs in relation to use of the device must be reported to the manufacturer and a competent authority of the Member State in which the user and / or patient resides.

DECLARATION OF CONFORMITY:

We are solely responsible for declaring that the Medical Devices mentioned in this statement are of Low-Risk Class (Class I) and comply with the requirements of the European Regulation 2017/745 and where appropriate, the standards and legislation referred to.

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MAINTENANCE, CLEANING & DISINFECTION INFORMATION:

To guarantee a safe use and adequate standard of hygiene, the user should perform these procedures before every use.
The user must make certain of the structural integrity of the device and its components. For cleaning and sanitation procedures, follow the steps below:

1. **Gloves should be worn.**
2. **Remove the evident dirt before to proceed with the most accurate cleaning.**
3. **Use water or mild detergent for the cleaning procedure.**
4. **Wipe with a dry and clean cloth.**
5. **DO NOT leave the product wet.**
6. **Do not use chemical cleaners to clean the frame. This could cause damage to the surface of the product.**
7. **If you need to disinfect the product, use a common, mild disinfectant.**
8. **Protect the item from scratches, cuts and punctures.**

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UNFOLDING THE WHEELCHAIR / AUSKLAPPEN DES ROLLSTUHLS

Gently hold the seat frame with one hand and the backrest frame with the other. Afterwards open with both arms outwards simultaneously and fasten the safety buckle downwards. Finaly turn the footrest and armrests downwards. In order to close follow the reverse procedure.

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
PARKING BRAKESFESTSTELLBREMSEN


Parking brakes (see picture) are used to completely immobilize the wheelchair and prevent it from rolling. To engage the parking brake, apply upward pressure while pressing the attendant’s brakes, and then release the drive brakes once you hear the distinctive click. To disengage the parking brakes, follow the reverse procedure.

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Never engage the parking brakes while you are moving to slow down the wheelchair.
Keep your fingers away from movable parts when using the brake, always keep your hand on the brake lever.

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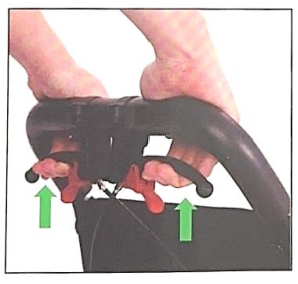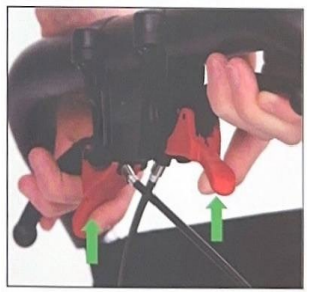PULL HANDLE/ ZUGGRIFF

Press the button on the special pull handle located under the seat and pull out to use it as a travel transfer wheelchair or easily place it in the trunk of the car. It also has a special pocket on the back of the backrest. It can also be easily transported on the plane.

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WARRANTY TERMS AND CONDITIONS:

The product is guaranteed for a period of two (2) years from the date of purchase.

The warranty covers manufacturing defects only. It does not cover defects or damage resulting from misuse, improper maintenance, modification, excessive use, or failure to comply with the product’s instructions for use. The warranty also does not cover parts that become worn through normal use.

In addition, the warranty does not cover damage or defects caused by natural disasters, unauthorized maintenance or repairs, power supply problems (where applicable), or transportation companies.

Any shipping or transportation costs related to a warranty claim shall be paid by the customer or dealer.

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Description of symbols that have been printed on the label, master carton and or user manual

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